

Ideal for basic counting and weighing applications, the LBK provides simple operation and a variety of modes. The LBK is easy to use, as it contains just five keys that are colour-coded to highlight frequently used functions. Sleek and compact design makes the LBK ideal for tight work spaces. Solidly constructed with durable ABS plastic housing, the LBK is equipped with a sealed keypad and large stainless steel pan that is easy to remove for cleaning.

The LBK provides the best value for quick, basic weighing tasks



Features

- Vivid, backlit LCD easily visible in any lighting conditions
- Colour-coded keys facilitate quick recognition of the most frequently used buttons
- Level indicator and adjustable feet ensure proper balance setup for optimum weighing results
- Large, grade 304 stainless steel pan allows easy cleaning
- Lightweight design simplifies transport and portability
- Solid ABS housing is durable for industrial use
- Audible overload alarm warns when capacity is exceeded
- External calibration allows for verification and adjustment with weights
- Parts counting with preset sample sizes
- Zero-tracking feature ensures display returns to zero reading
- AC adapter included
- Rechargeable battery included for operation almost anywhere
- Programmable auto power-off to save energy

Applications

- Weighing
- Parts counting
- Percentage weighing

About ADAM:


For 50 years, professionals worldwide have relied on Adam Equipment for an extensive selection of dependable and affordable weighing equipment. Adam provides the right balance of speed, performance and value.

Model	LBK 6
Capacity	6kg
Readability	1g
Repeatability (S.D.)	1g
Linearity (+/-)	2g
Pan Size	250x180mm
Weighing Units	g, kg, lb, lb:oz, oz
Stabilisation Time (sec)	2
Calibration	External calibration
Display	Backlit LCD with 24mm -high digits
Power Supply	Mains power cord / Internal rechargeable battery
Operating Temperature	0° to 40°C
Housing	ABS plastic
Overall Dim.	250x307x110mm (wxdxh)
Net Weight	1.8kg

Accessories

Item number	Description
700200056	In-use cover
700200063	In-use cover (pack of 5)
700200064	In-use cover (pack of 10)
700660290	Calibration certificate

Key Accessories

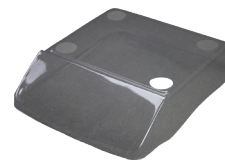

In-use cover
In-use covers provide an easy and affordable way to help keep the balance clean and protected.

Supplied by